

# Evaluation of Congo Red Staining Kit to Determine Proteinuria in Preeclampsia

JL RODRÍGUEZ CHÁVEZ<sup>1</sup>, EK FUENTES GUTIÉRREZ<sup>2</sup>, MJ ANGELESVÁZQUEZ<sup>3</sup>, H MENDIETA ZERÓN<sup>4</sup>

## ABSTRACT

**Introduction:** Preeclampsia is characterised by the increased excretion of misfolding protein, attracted by a Diazo dye denominated congo red.

**Aim:** To compare congo red staining kit versus a urocolor test to determine proteinuria in patients with preeclampsia.

**Materials and Methods:** This was a prospective, transversal study, performed between August 2010 to July 2011, including pregnant and puerperal women with the clinical diagnosis of preeclampsia. From all patients, a single sample of 24 hours urine was tested with congo red to determine the presence of proteinuria, comparing the qualitative results of the urine

dipstick for proteinuria and the gold standard determination of protein in 24 hours urine collection. Differences among nominal conditions were evaluated by the Mann-Whitney U-test.

**Results:** A total of 50 pregnant and puerperal women were included in this study. The statistical comparison between 2 vs. 3, 2 vs. 4 and 3 vs. 4, congo red positive fields reached significant difference ( $p \le 0.001$ ) in relation to the mean proteinuria values for each one. In contrast, with the urocolor test there were no differences per number of positive crosses.

**Conclusion:** The congo red has a high selectivity for betafolded proteins in preeclampsia being a promising screening test of high value to the clinician.

## INTRODUCTION

During pregnancy, protein excretion is increased due to a decrease in the proximal tubule reabsorption of Low-Molecular-Weight (LMW) proteins filtered by the glomerulus. Importantly, even healthy pregnant women can register proteinuria of >300 mg/dL [1].

Quantification of proteinuria in a 24 hours urine sample is mandatory in the comprehensive assessment of women with hypertensive disorders during pregnancy and, in some cases, it is necessary to establish the diagnosis of preeclampsia and its severity [2,3].

Significant proteinuria is defined by the National High Blood Pressure Education Program (NHBPEP) working group as an excretion of >300 mg of protein in a 24 hours urine sample. The same working group also states that proteinuria of 2 gm/day or greater involves increased disease severity, and this figure is set as a criterion for defining severe preeclampsia [4]. In contrast, the American College of Obstetrics and Gynaecology (ACOG) establishes a proteinuria of 5 gm/day or more for the diagnosis of severe preeclampsia [5].

Quantification of protein in a 24 hours urine collection is considered the gold standard for determining the presence of significant proteinuria in pregnancy [6]. However, the collection of urine every 24 hours is uncomfortable for the patient, time-consuming, and is subject to errors due to inadequate collection, delays the diagnosis for at least 24 hours, and it has been reported that fewer than one half of the sick women hospitalised with diagnosis of preeclampsia have a 24 hours urine collection. Due to the disadvantage of collecting 24 hours urine, other options have been proposed for the diagnosis of significant proteinuria in pregnancy. These other tests are fast and reliable, for example, the use of labstix or a test strip, measurement of proteinuria in short periods of urine collection, and the protein/urine creatinine ratio [7-9].

Protein three-dimensional conformation is the result of unfolded polypeptides that goes through a molecular folding process for which several models have been proposed. The misfolded

Keywords: Congophilia, Correlation, Misfolded proteins, Pregnancy

proteins produce amorphous aggregation, a phenomenon determined by the melting temperature (Tm) [10]. Some protein misfolding diseases include Alzheimer's, Parkinson's and some neurodegenerative diseases [11]. Interestingly, preeclampsia is characterised by the increased excretion of misfolding protein, attracted by a Diazo dye enominated congo red [12]. By using an unbiased mass spectrometry proteomic profiling approach, the group of Buhimschi IA et al., discovered that urine of women with severe preeclampsia displayed a characteristic set of proteomic biomarkers consisting of nonrandom proteoforms of SERPINA1 and albumin [13]. SERPINA1 fragments have the propensity to misfold and aggregate into supramolecular structures [14], thus showing congophilia, a characteristic of protein misfolding [15]. The aim of this study was to compare congo red staining kit versus a urocolor test to determine proteinuria in patients with preeclampsia, looking for a more specific option with good correlation to the 24 hours urine collection test.

#### MATERIALS AND METHODS

This was a prospective, descriptive, cross-sectional, pilot study conducted at the "Mónica Pretelini Sáenz" Maternal-Perinatal Hospital (HMPMPS), Health Institute of the State of Mexico (ISEM), Toluca from August 2010 to July 2011. The sample was set with a convenience non-probability sampling, based on the congo red tests available. Fifty pregnant women with >20 and <41 weeks gestation or in immediate puerperium exhibiting the clinical criteria of preeclampsia were enrolled in the study. The diagnosis of preeclampsia was consistent with the recommendations of the ACOG [5]. Women with a history of preeclampsia, Urinary Tract Infection (UTI), kidney disease, or gestational hypertension were excluded from the study. Patients with incomplete medical files were discarded from the final analysis. The present study was approved by the Ethics Committee of the HMPMPS (Code: El/136) and was conducted in accordance with the Helsinki Declaration (Fortaleza,

Brazil). Written informed consent was obtained from all patients who participated in the study.

#### **General Information**

The primary source of information was collected through the review of medical records and taking into account the following: age; educational level; factors associated with a gestational hypertensive disorder; number of sexual partners; prenatal care; type of delivery of the infant; main indication for cesarean section, and maternal and perinatal outcomes.


#### **Congo Red**

The Highman's method was adjusted using the semiquantitative Congo Red Kit (catalogue #101641; Merck Millipore, USA) [16]. The Congo red staining principle is based on the formation of hydrogen bridge bonds with the carbohydrate component of the substrate. Congo red is an anionic dye and is capable of depositing itself in amyloid fibrils, which then exhibit a conspicuous dichroism under polarised light. The tissue stained with Congo red appears orangered under the transmitted-light microscope; under polarised light; however, the amyloid deposits show up as brilliant green doublerefraction images against a dark background. The staining process was conducted by the Clinical Laboratory Service, taking an average of 12 minutes. The reading was performed by two chemical pharmacobiologists and an expert in anatomopathology. Each slide was divided in four visual fields at the moment of the microscopic visualisation to give a report based on the positive quadrants (fields).

#### **Process**

In Highman's method, the staining was done with an alcoholic solution (in present case 2% isopropyl alcohol was used), followed by a differentiation with an alkaline solution (KOH). In the present process the urinary sediment was used as a specimen to study and therefore steps of the process were minimised.

The sample was collected after a centrifugation for 3 minutes, and then it was fixed with heat. The next step was staining in congo red solution, followed by rinsing and drying. The interpretation was estimating its presence and the distribution in fields (quadrants) of vision in the microscope [Table/Fig-1]. It was reported as a positive result to dark blue nuclei whose center (amyloid) was pink-red pigmented and the periphery metachromasia green when using polarised light. Simultaneously, the 24 hours urine collection, and two measurements with urocolor 11 test strips within the first 6 hours, using a chromatic scale were also performed.

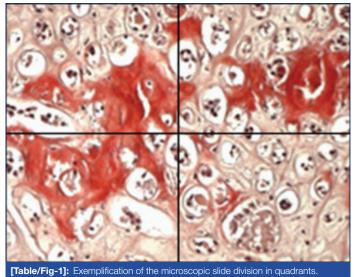

#### STATISTICAL ANALYSIS

Quantitative variables were expressed as mean±Standard Deviation (SD) and qualitative variables in frequency and percentage. The normality hypothesis was tested using the Kolmogorov-Smirnov test. Differences among nominal conditions were evaluated by the Mann-Whitney U-test. The intra and interobserver Kappa index were calculated for the identification of positive congophilia. In all cases, p≤0.05 was considered statistically significant. All tests were performed with the SPSS version 23.0 statistical software program (IBM SPPS 23, USA).

## RESULTS

#### **General Characteristics**

Considering 436 hospital admissions (0.88 admissions/day) into the Obstetrics Intensive Care Unit (OICU) during the study frame time, preeclampsia represented a total of 223 cases (51.14%) (0.44 admissions/day). From all this population only 50 women were selected. Mean age was 27.4±7.2 years, with an average hospital stay of 3.6±1.5 days and with three organ dysfunctions. The [Table/ Fig-2] shows the general characteristics of the studied population.

Variable	Value					
Age (years)	27.4±7.2					
Gestational age (weeks)	33.96±2.5					
Mean arterial pressure (mmHg)	111.6±24.3					
Proteinuria (mg/dL)	3076.2±1677.6					
Stay in the obstetrical ICU (days)	3.6±1.5					
Table/Fig-2]: General characteristics of the patients.						

In relation to the obstetrical data, 1 (28%) were primiparous and the average number of pregnancies was 2.72. Women were found at an obstetric gestational age of 33.96 weeks (early preeclampsia) and 45 (90%) had risk factors. The [Table/Fig-3] shows the distribution by age. Analysis of organ dysfunction was based on the Sequential Organ Failure Assessment (SOFA) classification and the main organ dysfunctions are described in [Table/Fig-4] [17]. At the moment of hospital discharge, >96% of patients had experienced remission of all dysfunctions. 26 (52%) of all women had fetal compromise. Fetal outcomes are described in [Table/Fig-5].

Age range (Years)	n (Percentage)					
15-20	9 (18%)					
21-25	12 (24%)					
26-30	14 (28%)					
31-35	6 (12%)					
36-40	7 (14%)					
≥40 years	2 (4%)					
Total	50 (100%)					
[Table/Fig.3]. Distribution of nationts per age range						

Organ failure	n (Percentage)				
Renal	50 (100%)				
Haematological	31 (62%)				
Hepatic	29 (58%)				
Cardiovascular	28 (56%)				
Neurological	13 (26%)				
Respiratory dysfunction	7 (14%)				
[Table/Fig-4]: Main organ dysfunctions.					

Outcome	Unit
APGAR score at 5 minutes (mean±SD)	7.4±1.8
Fetal weight (gm) (mean±SD)	1951.9±742.5
Intrauterine growth restriction (cases)	12 (24%)
Cesarean section (cases)	44 (88%)
<b>[Table/Fig-5]:</b> Fetal outcomes. SD: Standard deviation	

Congo red	N	Mean	Mean Standard deviation	95% Confidence interv Error for the mean		Standard Frror for the mea		Minimum	Maximum
fields			deviation		Lower limit	Upper limit			
2	12	1366.4	375.0	108.2	1128.0	1604.7	934.2	1892.0	
3	17	2664.5	471.0	114.2	2422.3	2906.7	1490.0	3349.7	
4	21	4386.7	1714.3	374.0	3606.3	5167.0	2763.1	9967.1	
Total	50	3076.2	1677.6	237.2	2599.5	3553.0	934.2	9967.1	
[Table/Fig-6]: Pro	Table/Fig-61: Proteinuria values per number of congo red positive fields.								

**Congo Red** 

Congophilia was present in two fields in 12 (24%), in three fields 34% (n=17), and in four fields, in 21 (42%). The intraobserver Kappa index was 0.87 (expected, 0.5) and the interobserver Kappa index was 0.740 {Standard Error (SE): 0.067; 95% Confidence Interval (95% CI), 0608-0872}.

In the whole population mean proteinuria value was 3,076.2 mg. In 4 (8%) protein determination was <1,000 mg, while in 10 (20%) this fell within the range of 1,001-2,000, 12 (24%) within the range of 2,001-3,000 mg, 15 (30%) within the nephrotic range between 3,001 and 4,000 mg, 4 (8%) exhibited between 4,001 and 5,000 mg and over 5,000 mg, 5 (10%). Visualising this information. The [Table/Fig-6] shows the mean proteinuria values per number of positive congo red fields and displays the information for urocolor 11 test strips [Table/Fig-7]. Contrastingly, although that 100% of cases were positive for congo red test none was 4+ positive in the urocolor 11 strip despite signs of nephrotic-range proteinuria [Table/Fig-8].

Kolmogorov-Smirnov normality test was used and evidenced a non-Gaussian distribution for proteinuria per group of red congo positive identification. In view of the results, Mann-Whitney U-test was used to compare the proteinuria values [Table/Fig-9].

On the other hand, the mean proteinuria values per number of crosses in the urocolor strip are shown in [Table/Fig-10]. Immediately stands out the fact that the proteinuria did not increase in relation to the crosses and the comparison was not significant [Table/Fig-11].

Congo red	Uro					
positive fields	None	None One Two		Three	Total	
2	6	6	0	0	12	
3	9	5	2	1	17	
4	6	12	3	0	21	
Total	21	23	5	1	50	
[Table/Fig-]	71: Number of	urocolor positi	ve crosses pe	r conao red po	ositive fields.	

	Positive congo Urolabstix					Proteinuria in
Case	microscope slides	red fields in the microscope slides 1+ 2+ 3+ 4+		24 hours urine collection (mg)		
1	3	0	0	0	0	2865
2	2	0	0	0	0	1765
3	3	0	0	0	0	2716
4	2	0	0	0	0	958.4
5	2	1	0	0	0	1892
6	3	0	0	0	0	2789
7	4	0	0	0	0	3209
8	4	1	0	0	0	9967.1
9	2	0	0	0	0	1243.6
10	3	0	0	0	0	2345.7
11	4	1	0	0	0	3522.2
12	3	0	0	0	0	2891
13	2	1	0	0	0	934.2
14	4	0	0	0	0	6540

Journal of Clinical and Diagnostic Research. 2018 May, Vol-12(5): QM01-QM05
---

J

1531003340164010038191730100278218410003276119401003228.120410003298.121300102952.42220000185623410003248.524401002850.224401003248.525410002840.22630100332928210003345.330410003345.331310003345.332210003325334100033253431000332535400003325364100033253721000332536400003325372100033253821							
17       3       0       1       0       0       2782         18       4       1       0       0       0       3276         19       4       0       1       0       0       3226         19       4       0       1       0       0       3202.1         20       4       1       0       0       0       3298.1         21       3       0       0       0       0       3298.1         22       2       0       0       0       0       3298.1         23       4       1       0       0       4563.2         24       4       0       1       0       0       3248.5         25       4       1       0       0       3248.5         26       3       0       1       0       0       3349.2         27       4       1       0       0       1455.3         28       2       1       0       0       2314.6         30       4       1       0       0       2314.6         32       2       1       0       0       0 <td>15</td> <td>3</td> <td>1</td> <td>0</td> <td>0</td> <td>0</td> <td>3240</td>	15	3	1	0	0	0	3240
18         4         1         0         0         0         3276           19         4         0         1         0         0         3202.1           20         4         1         0         0         0         3202.1           20         4         1         0         0         0         3295.1           21         3         0         0         1         0         295.4           22         2         0         0         0         1856           23         4         0         1         0         0         3248.5           24         4         0         1         0         0         2763.1           25         4         1         0         0         0         3249.2           26         3         0         1         0         0         3345.3           29         3         1         0         0         3456.2           31         3         1         0         0         3345.3           32         2         1         0         0         3325           33         4         1	16	4	0	1	0	0	3819
19         4         0         1         0         0         3202.1           20         4         1         0         0         0         3298.1           21         3         0         0         1         0         295.4           22         2         0         0         0         0         1856           23         4         1         0         0         0         3248.5           24         4         0         1         0         0         2763.1           26         3         0         1         0         0         2840.2           27         4         1         0         0         3329           28         2         1         0         0         3345.3           30         4         1         0         0         3345.3           31         3         1         0         0         1356.2           32         2         1         0         0         3345.3           33         4         1         0         0         3325           34         3         1         0         0	17	3	0	1	0	0	2782
20         4         1         0         0         0         3298.1           21         3         0         0         1         0         2952.4           22         2         0         0         0         0         1856           23         4         1         0         0         0         3248.5           24         4         0         1         0         0         3248.5           25         4         1         0         0         2763.1           26         3         0         1         0         0         3240.2           27         4         1         0         0         3329         3329           28         2         1         0         0         33329         345.3           30         4         1         0         0         0         3345.3           31         3         1         0         0         0         3345.3           32         2         1         0         0         0         3345.3           33         4         1         0         0         0         3345.3 <tr< td=""><td>18</td><td>4</td><td>1</td><td>0</td><td>0</td><td>0</td><td>3276</td></tr<>	18	4	1	0	0	0	3276
21         3         0         0         1         0         2952.4           22         2         0         0         0         0         1856           23         4         1         0         0         0         32455           24         4         0         1         0         0         3248.5           25         4         1         0         0         2763.1           26         3         0         1         0         0         2840.2           27         4         1         0         0         3329           28         2         1         0         0         3345.3           30         4         1         0         0         3345.3           31         3         1         0         0         3345.3           32         2         1         0         0         2314.6           33         4         1         0         0         3325           34         3         1         0         0         3325           35         4         0         0         0         3459	19	4	0	1	0	0	3202.1
22         2         0         0         0         0         0         1856           23         4         1         0         0         0         463.2           24         4         0         1         0         0         3248.5           25         4         1         0         0         2763.1           26         3         0         1         0         0         2840.2           27         4         1         0         0         3329           28         2         1         0         0         3345.3           30         4         1         0         0         3456.2           31         3         1         0         0         1356           32         2         1         0         0         1356           33         4         1         0         0         3325           34         3         1         0         0         3325           37         2         1         0         0         3459           40         0         0         0         0         3459 <td< td=""><td>20</td><td>4</td><td>1</td><td>0</td><td>0</td><td>0</td><td>3298.1</td></td<>	20	4	1	0	0	0	3298.1
23         4         1         0         0         0         44663.2           24         4         0         1         0         0         3248.5           25         4         1         0         0         0         2763.1           26         3         0         1         0         0         2840.2           27         4         1         0         0         0         3329           28         2         1         0         0         0         3345.3           30         4         1         0         0         0         3345.3           30         4         1         0         0         0         3345.3           31         3         1         0         0         0         3345.3           32         2         1         0         0         0         3345.3           33         4         1         0         0         0         3325           34         3         1         0         0         0         3325           37         2         1         0         0         34562	21	3	0	0	1	0	2952.4
24         4         0         1         0         0         3248.5           25         4         1         0         0         0         2763.1           26         3         0         1         0         0         2840.2           27         4         1         0         0         0         3329           28         2         1         0         0         0         3345.3           29         3         1         0         0         0         3345.3           30         4         1         0         0         0         3345.3           30         4         1         0         0         0         3345.3           31         3         1         0         0         0         3345.3           32         2         1         0         0         0         3345.3           33         4         1         0         0         0         3325           34         3         1         0         0         0         3325           36         4         1         0         0         3459	22	2	0	0	0	0	1856
25         4         1         0         0         0         2763.1           26         3         0         1         0         0         2840.2           27         4         1         0         0         3329           28         2         1         0         0         0         1455.3           29         3         1         0         0         0         3345.3           30         4         1         0         0         0         4356.2           31         3         1         0         0         0         2314.6           32         2         1         0         0         0         5567           34         3         1         0         0         0         2344           35         4         0         0         0         3325           36         4         1         0         0         980           38         2         1         0         0         3459           40         4         0         0         0         3459           41         4         0         0         0<	23	4	1	0	0	0	4563.2
26         3         0         1         0         0         2840.2           27         4         1         0         0         0         3329           28         2         1         0         0         0         3329           29         3         1         0         0         0         3345.3           30         4         1         0         0         0         4356.2           31         3         1         0         0         0         2314.6           32         2         1         0         0         0         2314.6           32         2         1         0         0         0         2314.6           34         3         1         0         0         0         2344           35         4         0         0         0         3325         333           36         4         1         0         0         0         34562           36         4         1         0         0         0         3459           40         4         0         0         0         3459 <td< td=""><td>24</td><td>4</td><td>0</td><td>1</td><td>0</td><td>0</td><td>3248.5</td></td<>	24	4	0	1	0	0	3248.5
27         4         1         0         0         0         3329           28         2         1         0         0         0         1455.3           29         3         1         0         0         0         3345.3           30         4         1         0         0         0         4356.2           31         3         1         0         0         0         2314.6           32         2         1         0         0         0         6577           34         3         1         0         0         0         2344           35         4         0         0         0         3325           37         2         1         0         0         980           38         2         1         0         0         1778           39         4         1         0         0         3459           40         4         0         0         0         3459           41         4         0         0         0         3298           42         3         1         0         0         3298 </td <td>25</td> <td>4</td> <td>1</td> <td>0</td> <td>0</td> <td>0</td> <td>2763.1</td>	25	4	1	0	0	0	2763.1
28         2         1         0         0         0         1455.3           29         3         1         0         0         0         3345.3           30         4         1         0         0         0         4356.2           31         3         1         0         0         0         2314.6           32         2         1         0         0         0         6577           34         3         1         0         0         0         2344           35         4         0         0         0         4562           36         4         1         0         0         980           37         2         1         0         0         980           38         2         1         0         0         980           38         2         1         0         0         1778           39         4         1         0         0         3459           40         4         0         0         0         3459           41         4         0         0         0         3298 <tr< td=""><td>26</td><td>3</td><td>0</td><td>1</td><td>0</td><td>0</td><td>2840.2</td></tr<>	26	3	0	1	0	0	2840.2
29         3         1         0         0         0         3345.3           30         4         1         0         0         0         4356.2           31         3         1         0         0         0         2314.6           32         2         1         0         0         0         1356           33         4         1         0         0         0         6577           34         3         1         0         0         0         2344           35         4         0         0         0         4562           36         4         1         0         0         980           38         2         1         0         0         980           38         2         1         0         0         1778           39         4         1         0         0         3459           40         4         0         0         0         2347.9           43         4         0         0         0         3298           44         4         0         0         0         1490 <t< td=""><td>27</td><td>4</td><td>1</td><td>0</td><td>0</td><td>0</td><td>3329</td></t<>	27	4	1	0	0	0	3329
30         4         1         0         0         0         4356.2           31         3         1         0         0         0         2314.6           32         2         1         0         0         0         1356           33         4         1         0         0         0         2314.6           32         2         1         0         0         0         1356           33         4         1         0         0         0         2314.6           34         3         1         0         0         0         2344           35         4         0         0         0         3325           36         4         1         0         0         980           38         2         1         0         0         1778           39         4         1         0         0         3459           40         4         0         0         0         2347.9           41         4         0         0         0         3298           44         0         0         0         0         1490	28	2	1	0	0	0	1455.3
31         3         1         0         0         0         2314.6           32         2         1         0         0         0         1356           33         4         1         0         0         0         6577           34         3         1         0         0         0         2344           35         4         0         0         0         0         4562           36         4         1         0         0         0         3325           37         2         1         0         0         0         980           38         2         1         0         0         0         1778           39         4         1         0         0         0         3459           40         4         0         0         0         3459           41         4         0         0         0         32347.9           43         4         0         0         0         3298           44         4         0         0         0         3298           44         0         0         0	29	3	1	0	0	0	3345.3
32         2         1         0         0         0         1356           33         4         1         0         0         0         6577           34         3         1         0         0         0         2344           35         4         0         0         0         0         4562           36         4         1         0         0         0         3325           37         2         1         0         0         0         980           38         2         1         0         0         0         3459           40         4         1         0         0         0         3459           41         4         0         0         0         3459           42         3         1         0         0         3287           43         4         0         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         1490           46         3         0         0         0	30	4	1	0	0	0	4356.2
33         4         1         0         0         0         6577           34         3         1         0         0         0         2344           35         4         0         0         0         0         2344           35         4         0         0         0         0         4562           36         4         1         0         0         0         3325           37         2         1         0         0         0         980           38         2         1         0         0         0         3459           40         4         0         0         0         3459           40         4         0         0         0         3459           41         4         0         0         0         3456           41         4         0         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         1490           46         3         0         0         0         21011	31	3	1	0	0	0	2314.6
34         3         1         0         0         0         2344           35         4         0         0         0         0         4562           36         4         1         0         0         0         3325           37         2         1         0         0         0         980           38         2         1         0         0         0         3459           40         4         0         0         0         3459           40         4         0         0         0         3257           40         4         0         0         0         3459           40         4         0         0         0         3456           41         4         0         0         0         2347.9           42         3         1         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         1490           46         3         0         0         0         210.1           48	32	2	1	0	0	0	1356
35         4         0         0         0         0         4562           36         4         1         0         0         0         3325           37         2         1         0         0         0         980           38         2         1         0         0         0         1778           39         4         1         0         0         0         3459           40         4         0         0         0         4356           41         4         0         0         0         2347.9           42         3         1         0         0         2347.9           43         4         0         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         1490           46         3         0         0         0         2210.1           47         2         0         0         0         3349.7           48         2         0         0         0         3349.7           5	33	4	1	0	0	0	6577
36         4         1         0         0         0         3325           37         2         1         0         0         0         980           38         2         1         0         0         0         1778           39         4         1         0         0         0         3459           40         4         1         0         0         0         4356           41         4         0         0         0         2347.9           42         3         1         0         0         2347.9           43         4         0         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         1490           46         3         0         0         0         1490           47         2         0         0         0         3349.7           50 <td>34</td> <td>3</td> <td>1</td> <td>0</td> <td>0</td> <td>0</td> <td>2344</td>	34	3	1	0	0	0	2344
37         2         1         0         0         0         980           38         2         1         0         0         0         1778           39         4         1         0         0         0         3459           40         4         1         0         0         0         4356           41         4         0         0         0         2347.9           42         3         1         0         0         0         2347.9           43         4         0         0         0         3298         3298           44         4         0         0         0         0         3298           44         4         0         0         0         3298         33           44         4         0         0         0         3298         33         34           46         3         0         0         0         1490         34           47         2         0         0         0         349.7         349         33         0         0         0         3349.7           50         3 <t< td=""><td>35</td><td>4</td><td>0</td><td>0</td><td>0</td><td>0</td><td>4562</td></t<>	35	4	0	0	0	0	4562
38         2         1         0         0         0         1778           39         4         1         0         0         0         3459           40         4         1         0         0         0         3459           41         4         0         0         0         0         4356           41         4         0         0         0         0         5890           42         3         1         0         0         0         2347.9           43         4         0         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         3298           44         4         0         0         0         3298           45         3         0         0         0         1490           46         3         0         0         0         1490           48         2         0         0         0         3349.7           50         3         0         0         0         2474	36	4	1	0	0	0	3325
39         4         1         0         0         0         3459           40         4         1         0         0         0         3459           41         4         0         0         0         0         4356           41         4         0         0         0         0         5890           42         3         1         0         0         0         2347.9           43         4         0         0         0         0         3298           44         4         0         0         0         0         5560.2           45         3         0         0         0         0         1490           46         3         0         0         0         0         2210.1           47         2         0         0         0         0         198.2           48         2         0         0         0         0         3349.7           50         3         0         0         0         0         2474           Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours         12474         124 <t< td=""><td>37</td><td>2</td><td>1</td><td>0</td><td>0</td><td>0</td><td>980</td></t<>	37	2	1	0	0	0	980
40         4         1         0         0         0         4356           41         4         0         0         0         0         23690           42         3         1         0         0         0         2347.9           43         4         0         0         0         0         3298           44         4         0         0         0         0         3298           44         4         0         0         0         0         3298           44         4         0         0         0         0         3298           45         3         0         0         0         0         3298           46         3         0         0         0         0         1490           47         2         0         0         0         1198.1           48         2         0         0         0         3349.7           50         3         0         0         0         2474           [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours         349.7         349.7         349.7	38	2	1	0	0	0	1778
41         4         0         0         0         0         0         5890           42         3         1         0         0         0         2347.9           43         4         0         0         0         0         3298           44         4         0         0         0         0         3298           44         4         0         0         0         0         3298           45         3         0         0         0         0         3298           46         3         0         0         0         1490         1490           46         3         0         0         0         0         2210.1           47         2         0         0         0         0         1198.1           48         2         0         0         0         980.2           49         3         0         0         0         3349.7           50         3         0         0         0         2474           [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours         24         34         34	39	4	1	0	0	0	3459
42         3         1         0         0         0         2347.9           43         4         0         0         0         0         3298           44         4         0         0         0         0         3298           44         4         0         0         0         0         5560.2           45         3         0         0         0         0         1490           46         3         0         0         0         0         2210.1           47         2         0         0         0         198.1           48         2         0         0         0         980.2           49         3         0         0         0         2474           [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours         2474	40	4	1	0	0	0	4356
43       4       0       0       0       0       3298         44       4       0       0       0       0       5560.2         45       3       0       0       0       0       1490         46       3       0       0       0       0       2210.1         47       2       0       0       0       0       1198.1         48       2       0       0       0       0       980.2         49       3       0       0       0       0       3349.7         50       3       0       0       0       0       2474         [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours	41	4	0	0	0	0	5890
44         4         0         0         0         0         0         5560.2           45         3         0         0         0         0         0         1490           46         3         0         0         0         0         0         2210.1           47         2         0         0         0         0         198.1           48         2         0         0         0         0         980.2           49         3         0         0         0         0         3349.7           50         3         0         0         0         2474           [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours	42	3	1	0	0	0	2347.9
45         3         0         0         0         0         1490           46         3         0         0         0         0         1490           46         3         0         0         0         0         2210.1           47         2         0         0         0         0         1198.1           48         2         0         0         0         0         980.2           49         3         0         0         0         0         3349.7           50         3         0         0         0         2474           [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours	43	4	0	0	0	0	3298
46       3       0       0       0       0       2210.1         47       2       0       0       0       0       1198.1         48       2       0       0       0       0       980.2         49       3       0       0       0       0       3349.7         50       3       0       0       0       2474         [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours	44	4	0	0	0	0	5560.2
47       2       0       0       0       0       1198.1         48       2       0       0       0       0       980.2         49       3       0       0       0       0       3349.7         50       3       0       0       0       0       2474         [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours	45	3	0	0	0	0	1490
48         2         0         0         0         0         980.2           49         3         0         0         0         0         3349.7           50         3         0         0         0         0         2474           (Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours         24 hours         24 hours	46	3	0	0	0	0	2210.1
49         3         0         0         0         0         3349.7           50         3         0         0         0         0         2474           [Table/Fig-8]:         Congo red positive fields, urocolor test and proteinuria in 24 hours	47	2	0	0	0	0	1198.1
50         3         0         0         0         2474           [Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours	48	2	0	0	0	0	980.2
[Table/Fig-8]: Congo red positive fields, urocolor test and proteinuria in 24 hours	49	3	0	0	0	0	3349.7
	50	3	0	0	0	0	2474
	[Table/I	Fig-8]: Congo red posit					

and the second				
urine collection.	Irine	col	lect	ion.

Groups of comparison per congo red positive fields	Mann-Whitney U-test	W of Wilcoxon	z	Sig.				
2 vs. 3	4	82	-4.34	≤0.001				
2 vs. 4	0.0	78	-4.715	≤0.001				
3 vs. 4	27	180	-4.448	≤0.001				
[Table/Fig-9]: Comparison of proteinuria per congo red positive fields.								

[Table/Fig-9]: Comparison of proteinuria per congo red positive fiel Sig: Significance

Number of	N	N Mean Standard Error for the me		Standard Error		Mean Standard Error			Minimum	Maximum
crosses			deviation		Lower limit	Upper limit				
none	21	2866.2	1595.8	348.2	2139.8	3592.6	958.4	6540.0		
1	23	3251.2	1960.7	408.8	2403.3	4099.1	934.2	9967.1		
2	5	3178.3	414.5	185.3	2663.6	3693.0	2782.0	3819.0		
3	1	2952.4					2952.4	2952.4		
Total	50	3076.2	1677.6	237.2	2599.5	3553.0	934.2	9967.1		
Table/Fig-101. M	lean proteinuria valu	ies per number of o	rosses in the Lirocol	or strin						

[Table/Fig-10]: Mean proteinuria values per number of crosses in the Urocolor strip.

Groups of comparison per urocolor positive crosses	Mann-Whitney U-test	W of Wilcoxon	z	Sig.
None vs. 1	203	434	-0.905	0.366
None vs. 2	36	267	1.073	0.283
None vs. 3	7	238	-0.552	0.581
1 vs. 2	56	332	-0.090	0.928
1 vs. 3	10	11	-0.217	0.828
2 vs. 3	2	3	-0.293	0.770
[Table/Fig-11]: Comparison of proteinuria per urocolor positive crosses. Sig: Significance				

## DISCUSSION

As noted, the percentage of caesarean sections exceeded the expected international goals, with 44 cases (88%). This is a national trend that has been widely discussed [18,19].

Today, many hospitals include the use of urine dipstick in the comprehensive evaluation of women with hypertensive disorders of pregnancy because it comprises part of international criteria for preeclampsia [20,21]. The urine dipstick is inexpensive and easy to use, and the result of the application of this resource is that it determines an approximate protein concentration, but not an absolute value. The patient's hydration status can influence the result; in fact, a false-negative result may occur in urine samples with a density of <1.010, or in the presence of high salt concentrations, with a very acidic proteinuria pH or proteinuria not due to the presence of albumin. All of these factors render the urine dipstick poorly correlated with the 24 hours proteinuria determination, leading to a low diagnostic performance as an indicator of significant proteinuria [22]. The analysis of several studies that have investigated the reliability of the urine dipstick have found low sensitivity and specificity, resulting in a diagnosis of preeclampsia that is poorly sustained or underdiagnosed [23,24]. This last sentence is accordance with our results as in some severe cases of proteinuria, identified both by 24 hours urine collection and congo red, the urine dipstick was negative. Therefore, with the urocolor strip there was poor association with disease severity ranges.

Although, determination of proteinuria by dipstick is deficient compared with the proteinuria value in the 24 hours urine sample [25], it is widely used. Even more, the absence of proteinuria on dipstick testing may have ruled out the presence of significant proteinuria in patients. Another alternative for the presence-of-proteinuria feature has been collecting urine for short periods (8 or 12 hours), which shares limitations and problems with the collection of 24 hours urine [26]. In this current manuscript, 38% of the patients had proteinuria  $\geq$ 3 gm/dL, illustrating the severity of present cases.

In a pilot study based on the congophilia test in urine, Buhimschi IA et al., found that 4/35 (11%) of asymptomatic women were positive for congo red prior to the development of preeclampsia [12]. These findings indicate the possibility of developing a test to predict and diagnose preeclampsia that only requires a few drops of urine, according to the author.

A previous study revealed a Detection Rate (DR) of 93% and a

False-Positive Rate (FPR) of 4%. However, with first-trimester urine samples, DR was 33.3, 16.1, and 20% for early, late, and all cases of preeclampsia, respectively, with a 12.8% FPR. The Odds Ratio (OR) for preeclampsia was superior to Body Mass Index (BMI), and Mean Arterial Blood Pressure (MAP), but inferior to previous preeclampsia and black ethnicity. In the first trimester, the Congo Red urine test adds accuracy for preeclampsia prediction in obese, black women who had previous preeclampsia and above-average mean arterial pressure [27].

In the present pilot study, the positivity with congo red was 100%, but we have to be aware that the observed variability in the positive fields is an issue that can be used for quality and interobserver analysis.

Recently, Parikh L et al., in a study with 101 women over 14 weeks gestation undergoing evaluation for preeclampsia and diabetes proved that the Congo Red Retention (CRR) was not increased in diabetic gravidas without preeclampsia even in the presence of baseline nephropathy. They also determined that a CRR over 15% correctly classified all women with preeclampsia and adverse pregnancy outcomes [28].

It is clear from this study that there is an increase in proteinuria values, with the Congo Red staining from two to three positive fields in 1298.1 mg and 1722.2 mg from three to four positive fields. Furthermore, the statistical comparison between 2 vs. 3, 2 vs. 4 and 3 vs. 4 congo red positive fields at the microscope slide visualisation reached significant difference ( $p \le 0.001$ ) en relation to the mean proteinuria values for each one, contrasting with the lack of significance with the urocolor test, an expected result due to the heterogeneous proteinuria increase per number of positive crosses in the strip.

The cost of the congo red it is an important issue. Although, it is more expensive than the urocolor tests available in the market, as one can infer from the statistical analysis, the second option does not offer congruent data with the proteinuria values, thus putting patients at risk of a misdiagnosis and incorrect treatment. A handicap of the CRR is that this technique might not be able to differentiate between autoimmune diseases like lupus nephritis and a hypertensive disorder of pregnancy, and this condition applies for other kits available in the market [29-31].

Preeclampsia is still one of leading causes of maternal mortality worldwide; unfortunately, reliable and economically affordable biomarkers to predict it are elusive [32]. Herein, combining biomarkers may improve the prediction of pregnancy hypertension in the early stages of gestation and we believe CRR could add to the positive predictive value of a set of tests to diagnose preeclampsia at an early stage [33].

Despite the difficulty to explore the congophilia with the selected kit, the results lead us to think about the possibility of designing a kit or a test strip using congo red to determine proteinuria at low cost, using it at home and facilitate the early detection of a possible obstetric complication. In fact, Jonas SM et al., have developed a smart phone-based diagnostic for preeclampsia administering congo red [34]. Further research is required to explore the use of congophilia as a common diagnosis test of proteinuria in pregnancy.

#### LIMITATION

A clear limitation of this study is that it was not a randomised study; also, a safety matter of concern is that the congo red stain might be carcinogenic. Finally, a second step to contrast this test with a control group and different techniques to quantify proteinuria was not possible to perform as the used kit was discarded from the Mexican market.

## CONCLUSION

This study confirms that women with preeclampsia and chronic kidney disease without preeclampsia have elevated congophilia urine levels. In this regard, the CRR interpreted a positive fields in the slide which shows a direct relationship with the proteinuria levels. Thus, the congo red test point represents a new diagnostic tool for the clinician working in the obstetrics area, who are looking for a quick and highly reliable test. This is a test that can be performed in any context, even in third-world countries, and the results after 10-15 minutes are obtained. In summary, the overall assessment of the misfolded protein load by congo red staining is expected to be a simple, efficient and highly reproducible diagnostic test. The results are very encouraging to continue a new line of research. By contrast, the urocolor test is not reliable to make a decision for renal injury in pregnancy/puerperium.

### ACKNOWLEDGEMENTS

Authors thank the staff of the Clinical Laboratory of the HMPMPS for their help to manage the kit and Maggie Brunner M.A., for her excellent help with the English style correction.

#### REFERENCES

- Phillips JK, McBride CA, Hale SA, Solomon RJ, Badger GJ, Bernstein IM. Examination of prepregnancy and pregnancy urinary protein levels in healthy nulliparous women. Reprod Sci. 2017;24(3):407-12.
- [2] Morgan JL, Nelson DB, Roberts SW, Wells CE, McIntire DD, Cunningham FG. Association of baseline proteinuria and adverse outcomes in pregnant women with treated chronic hypertension. Obstet Gynecol. 2016;128(2):270-76.
- [3] Berry C, Atta MG. Hypertensive disorders in pregnancy. World J Nephrol. 2016;5(5):418-28.
- [4] Report of the National High Blood Pressure Education Program Working Group on High Blood Pressure in Pregnancy. Am J Obstet Gynecol. 2000;183(1):S1-22.
- [5] Roberts JM, August PA, Bakris G, Barton JR, Bernstein IM, Druzin M, et al. American College of Obstetricians and Gynecologists, Task Force on Hypertension in Pregnancy. Hypertension in pregnancy. Report of the American College of Obstetricians and Gynecologists' Task Force on Hypertension in Pregnancy. Obstet Gynecol. 2013;122(5):1122-31.
- [6] Henderson JT, Thompson JH, Burda BU, Cantor A, Beil T, Whitlock EP. Screening for preeclampsia: a systematic evidence review for the U.S. Preventive Services Task Force [Internet]. Rockville (MD): Agency for Healthcare Research and Quality (US); 2017 Apr.
- [7] Kieler H, Zettergren T, Svensson H, Dickman PW, Larsson A. Assessing urinary albumin excretion in pre-eclamptic women: which sample to use? BJOG. 2003;110(1):12-17.
- [8] Somanathan N, Farrell T, Galimberti A. A comparison between 24-h and 2-hour urine collection for the determination of proteinuria. J Obstet Gynaecol. 2003;23(4):378-80.
- Wikström AK, Wikström J, Larsson A, Olovsson M. Random albumin/creatinine ratio for quantification of proteinuria in manifest pre-eclampsia. BJOG. 2006;113(8):930-34.
- [10] Khan MV, Zakariya SM, Khan RH. Protein folding, misfolding and aggregation: A tale of constructive to destructive assembly. Int J Biol Macromol. 2018;112:217-29.

- JL Rodríguez Chávez et al., Congo Red Staining in Preeclampsia
- [11] Klaips CL, Jayaraj GG, Hartl FU. Pathways of cellular proteostasis in aging and disease. J Cell Biol. 2018;217(1):51-63.
- [12] Buhimschi IA, Nayeri UA, Zhao G, Shook LL, Pensalfini A, Funai EF, et al. Protein misfolding, congophilia, oligomerization, and defective amyloid processing in preeclampsia. Sci Transl Med. 2014;6(245):245ra92.
- [13] Buhimschi IA, Zhao G, Funai EF, Harris N, Sasson IE, Bernstein IM, et al. Proteomic profiling of urine identifies specific fragments of SERPINA1 and albumin as biomarkers of preeclampsia. Am J Obstet Gynecol. 2008;199(5):551. e1-16.
- [14] Huntington JA. Serpin structure, function and dysfunction. J Thromb Haemost. 2011;9 Suppl 1:26-34.
- [15] Elghetany MT, Saleem A. Methods for staining amyloid in tissues: A review. Stain Technol. 1988;63(4):201-12.
- [16] Available at: http://www.merckmillipore.com/MX/es/product/Congo-Redstaining-kit,MDA\_CHEM-101641?ReferrerURL=https%3A%2F%2Fwww. google.com.mx%2F. Date of access: 01/02/2018.
- [17] Jain S, Guleria K, Suneja A, Vaid NB, Ahuja S. Use of the Sequential Organ Failure Assessment score for evaluating outcome among obstetric patients admitted to the intensive care unit. Int J Gynaecol Obstet. 2016;132(3):332-36.
- [18] Estrada-Aguilar CG, Hernández-Carranza R, Gómez-Alonso C, Negrete-Paz M de L. [Incidence of cesareansection in a regional general hospital]. Rev Med Inst Mex Seguro Soc. 2012;50(5):517-22.
- [19] Martínez-Salazar GJ, Grimaldo-Valenzuela PM, Vázquez-Peña GG, Reyes-Segovia C, Torres-Luna G, Escudero-Lourdes GV. [Caesarean section: History, epidemiology, and ethics to diminish its incidence]. Rev Med Inst Mex Seguro Soc. 2015;53(5):608-15.
- [20] Baba Y, Yamada T, Obata-Yasuoka M, Yasuda S, Ohno Y, Kawabata K, et al. Urinary protein-to-creatinine ratio in pregnant women after dipstick testing: prospective observational study. BMC Pregnancy Childbirth. 2015;15:331.
- [21] Baba Y, Ohkuchi A, Usui R, Takahashi H, Matsubara S. Urinary protein-tocreatinine ratio indicative of significant proteinuria in normotensive pregnant women. J Obstet Gynaecol Res. 2016;42(7):784-88.
- [22] Amin SV, Illipilla S, Hebbar S, Rai L, Kumar P, Pai MV. Quantifying proteinuria in hypertensive disorders of pregnancy. Int J Hypertens. 2014;2014:941408.
- [23] Abebe J, Eigbefoh J, Isabu P, Okogbenin S, Eifediyi R, Okusanya B. Accuracy of urine dipsticks, 2-h and 12-h urine collections for protein measurement as compared with the 24-h collection. J Obstet Gynaecol. 2008;28(5):496-500.
- [24] Gangaram R, Naicker M, Moodley J. Accuracy of the spot urinary microalbumin:creatinine ratio and visual dipsticks in hypertensive pregnant women. Eur J ObstetGynecol Reprod Biol. 2009;144(2):146-48.
- [25] De Silva DA, Halstead AC, Côté AM, Sabr Y, von Dadelszen P, Magee LA. Urinary dipstick proteinuria testing: does automated strip analysis offer an advantage over visual testing? J Obstet Gynaecol Can. 2014;36(7):605-12.
- [26] Stout MJ, Conner SN, Colditz GA, Macones GA, Tuuli MG. The utility of 12-hour urine collection for the diagnosis of preeclampsia: A systematic review and metaanalysis. Obstet Gynecol. 2015;126(4):731-36.
- [27] Sammar M, Syngelaki A, Sharabi-Nov A, Nicolaides K, Meiri H. Can staining of damaged proteins in urine effectively predict preeclampsia? Fetal Diagn Ther. 2017;41(1):23-31.
- [28] Parikh L, Iqbal SN, Desale S, Buhimschi IA, Umans JG. 988: Urine congophilia in gravidas with diabetes and clinically diagnosed preeclampsia. Am J Obstet Gynecol. 2017;216(1).S556.
- [29] McCarthy FP, Adetoba A, Gill C, Bramham K, Bertolaccini M, Burton GJ, et al. Urinary congophilia in women with hypertensive disorders of pregnancy and preexisting proteinuria or hypertension. Am J Obstet Gynecol. 2016;215(4):464.e1-7.
- [30] Available at: https://www.sigmaaldrich.com/catalog/product/sigma/ ht60?lang=es&region=MX. Date of access: 04/02/2018.
- [31] Available at: http://www.abcam.com/congo-red-stain-kit-amyloid-stainab150663.html. Date of access: 04/02/2018.
- [32] Nissaisorakarn P, Sharif S, Jim B. Hypertension in pregnancy: defining blood pressure goals and the value of biomarkers for preeclampsia. Curr Cardiol Rep. 2016;18(12):131.
- [33] Sharma K, Singh R, Kumar M, Gupta U, Rohil V, Bhattacharjee J. First-trimester inflammatory markers for risk evaluation of pregnancy hypertension. J Obstet Gynaecol India. 2018;68(1):27-32.
- [34] Jonas SM, Deserno TM, Buhimschi CS, Makin J, Choma MA, Buhimschi IA. Smartphone-based diagnostic for preeclampsia: an mHealth solution for administering the Congo Red Dot (CRD) test in settings with limited resources. J Am Med Inform Assoc. 2016;23(1):166-73.

#### PARTICULARS OF CONTRIBUTORS:

- 1. PhD. Department of Intensive Care Unit, "Mónica Pretelini Sáenz" Maternal-Perinatal Hospital, Toluca, State of Mexico, Mexico.
- 2. M.D. Department of Intensive Care Unit, "Mónica Pretelini Sáenz" Maternal-Perinatal Hospital, Toluca, State of Mexico, Mexico.
- M.D. Department of Intensive Care Unit, "Mónica Pretelini Sáenz" Maternal-Perinatal Hospital, Toluca, State of Mexico, Mexico.
   PhD. Department of Intensive Care Unit, "Mónica Pretelini Sáenz" Maternal-Perinatal Hospital, Toluca, State of Mexico, Mexico.
- NAME, ADDRESS, E-MAIL ID OF THE CORRESPONDING AUTHOR:

#### Dr. H Mendieta Zerón,

Felipe Villanueva sur 1209, Toluca, State of Mexico, Mexico. E-mail: mezh\_74@yahoo.com

FINANCIAL OR OTHER COMPETING INTERESTS: None.

Date of Submission: Dec 07, 2017 Date of Peer Review: Jan 29, 2018 Date of Acceptance: Feb 21, 2018 Date of Publishing: May 01, 2018